

**VWR** 

Supplier Partnerships for Customer Solutions

**1.800.932.5000** **vwr.com**

Prices, product appearance and specifications are current at the time of printing, subject to change without notice. Availability for certain products may be limited by federal, state, provincial or local licensing requirements. All prices are in U.S. dollars unless otherwise noted. Offers valid in USA, void where prohibited by law or company policy, while supplies last. Visit [vwr.com](http://vwr.com) to view our privacy policy and additional disclaimers.

VWR, forms of VWR and the VWR logo and/or design are either registered trademarks®, trademarks™, or service marksSM of VWR International, Inc. in the United States and/or other countries. All other marks referenced herein are registered trademarks, trademarks or service marks of their respective owner(s). For a complete list of trademark owners, please visit [vwr.com](http://vwr.com)

©2007 VWR International, Inc. All rights reserved. Printed in U.S.A.

254503-001 Rev.A



## VWR symphony Meter Abridged User's Guide



**symphony™**

# Table of Contents

<b>Chapter I</b> Introduction	
Password Protected Methods	1
AUTO-READ™	1
Stirrer Control	2
<b>Chapter II</b> Display	
General Description	3
<b>Chapter III</b> Keypad	
General Description	5
Convenience Features	6
Audible Signals & Alarms	6
Automatic Shut-off	6
Icon Definitions	7
<b>Chapter IV</b> Preparation	
Installing the Power Adapter	9
Installing Batteries	10
Connecting the Electrodes	12
Turning on the Instrument	14
<b>Chapter V</b> Menu Setup	
Navigation Tips	15
General Menu Settings	16
Time and Date Setup	18
Continuous, Timed, or AUTO-READ™ Measurement Selection	20
Selecting Measurement Mode	22
Setting Calibration	23

## Chapter VI pH Technique

pH Setup .....	25
pH Calibration .....	27
pH Measurement .....	28

## Chapter VII Dissolved Oxygen Technique

DO Setup .....	29
DO Calibration .....	31
DO Measurement .....	33

## Chapter VIII Conductivity Technique

Conductivity Setup .....	35
Conductivity Calibration .....	37
Conductivity Measurement .....	39

## Chapter IX ISE Technique

ISE Setup Options .....	41
Preparation of Standards .....	43
ISE Calibration .....	44
ISE Measurement .....	45

## Chapter X Declaration of Conformity

WEEE Compliance .....	47
-----------------------	----

# Chapter I Introduction

This user guide has been abridged to include the most pertinent details. To download the extended user guide, visit our website at [www.vwr.com](http://www.vwr.com).

Congratulations! You have selected a VWR symphony meter designed for electrochemistry measurement in the field or in the laboratory.

- Single parameter meters provide single measurement of pH, dissolved oxygen, or conductivity.
- Dual parameter meters provide dual measurements of pH/dissolved oxygen, pH/conductivity, or pH/ISE (ion selective electrode).
- Multiple parameter meters provide multiple measurements of pH, ISE, dissolved oxygen, and conductivity.

Built to meet the demands of busy, multi-user laboratories or plant environments, all meters are microprocessor controlled aiding in the delivery of accurate and precise measurements. The waterproof portable meters can even withstand submersion for short periods of time without any negative effect on operation. To better meet the needs of users in environmental protection and control, food and beverage, pharmaceutical, and consumer product laboratories, the symphony meters include these key features:

- **Password Protected Methods** – Meter's memory will save up to ten custom measurements and calibrations for future reference. Password protection of each method eliminates any tampering with methods as multiple users access only the procedure most appropriate to their work.
- **AUTO-READ™** – The instrument starts a measurement and automatically prints or logs data when the reading becomes stable.

## Introduction


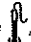
- **Stirrer Control** – Benchtop meters have a stirrer control for the Stirrer Probe and the AUTO-STIR™ BOD probe, eliminating the need for additional stir plates and stir bars.

An easy-to-use reference guide, attached to each meter, supports daily use.


Please read this manual thoroughly before using your benchtop or portable meter. Any use outside of these instructions could invalidate your warranty and/or cause permanent damage to the meter.

## Chapter II Display
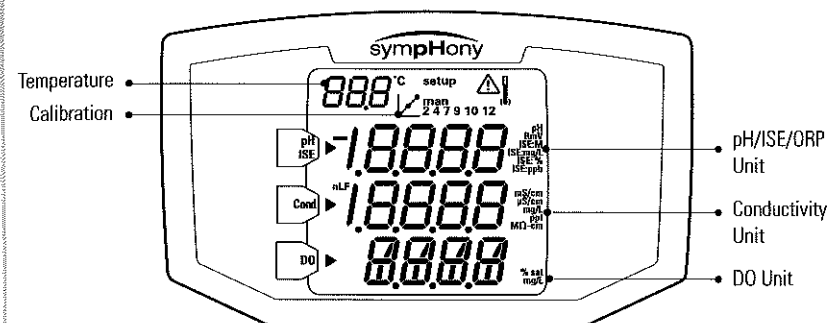
### General Description

Throughout a given process, the LCD on any sympHony meter provides **Temperature** and **Calibration** data. The **setup** only appears when the meter is in setup mode. The  indicates an error condition; when displayed with the , a sensor quality issue exists.


**AR** will be discussed in greater detail in **Chapter V, Menu Setup**.



### Multi-Parameter Meter

Depicted here is the primary LCD of the sympHony meter capable of multi-parameter measurement.

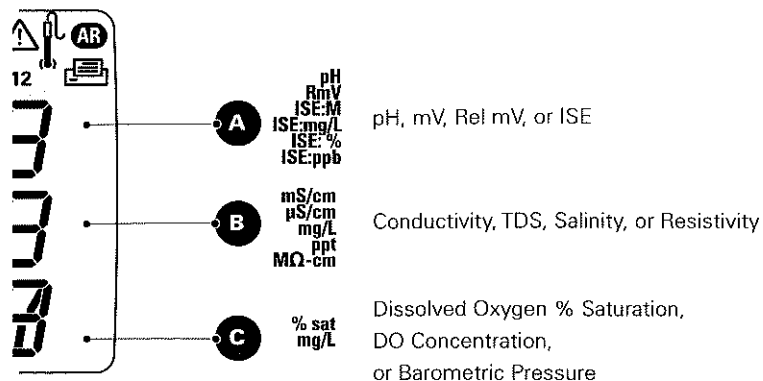
**Note:** the lowest 3 lines of data correspond to what is being measured. ▲




Units of measure, displayed on the right side of the screen, will flash until the reading is stable.

## Chapter III Keypad

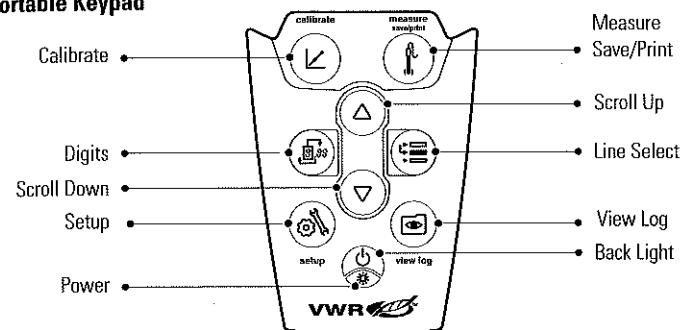
### General Description



Keypad layout is the same for all symphony meters. The portable meters have 9 keys. The benchtop meters have 10 keys due to the addition of a  button.

#### Portable Keypad



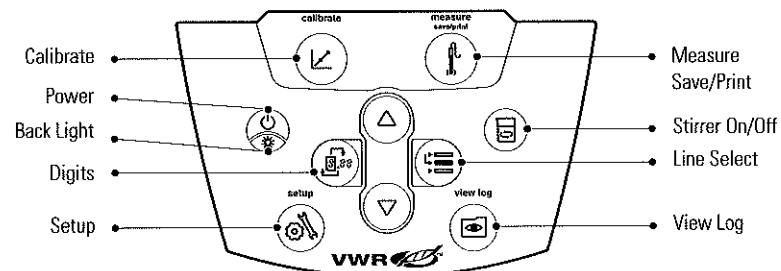


#### Benchtop Keypad


















The center area of the ergonomically designed keypad is the focal point for setting up and navigating the meter. In particular,  and  are used frequently to change LCD views.

#### Convenience features include:

- **Display Backlight** – A quick press of  will turn the backlight on and off. When the meter is operating on battery power, the backlight will automatically turn off after two minutes to conserve power. When batteries are low, the backlight will no longer turn on.
- **Automatic Shut-off** – All symphony meters will shut down after 20 minutes without a keypress. This maximizes battery power on portable meters, or benchtop meters being run on battery power.
- **Audible Signals** – The meter will beep whenever a key is pressed providing immediate verification that the user's input was received.
- **Visual Alarm Signals** – Flashing  and  indicate calibration settings need adjustment. For more detail, refer to subsequent sections of this manual discussing specific measurement techniques.

To become familiar with each function, please review these icon definitions.

## Icon Definitions

Key	Description	Key	Description
	<ul style="list-style-type: none"> <li>• Turns meter on, if shut off.</li> <li>• Toggles backlight on and off, if meter is on.</li> <li>• If meter is on, holding down  will turn off the meter.</li> </ul>		<ul style="list-style-type: none"> <li>• Changes measurement mode of the selected line.</li> <li>• Changes selected line in setup, methods and log view modes.</li> </ul>
			<ul style="list-style-type: none"> <li>• Edits the value of the flashing digit for setup, password entry, and calibration modes.</li> </ul>
	Switches arrow on left of the screen between 3 display lines to select and edit.		Changes selected digit to edit and moves decimal point when changing values in setup, password entry, and calibration.
	<ul style="list-style-type: none"> <li>• Starts calibration for currently selected line in measurement mode.</li> <li>• If arrow is pointed to the top line and the current units are pH,  will start a pH calibration.</li> <li>• Each time  is pressed in the calibration it will accept the current point value and move to the next cal point until the max # of cal pts are performed, then it will return to measurement mode.</li> </ul>		<ul style="list-style-type: none"> <li>• Prints and logs a measurement in continuous or timed measurement modes.</li> <li>• Prints, logs, and freezes screen when the reading becomes stable in AUTO-READ™ mode.</li> <li>• Exits setup menu and returns to measurement mode.</li> <li>• Accepts calibration and returns to measurement mode.</li> </ul>
	<ul style="list-style-type: none"> <li>• Enters setup menu starting with selected line in measurement mode:</li> <li>• If arrow is pointed to top line and current units are ISE, then  will enter the ISE setup screen.</li> </ul>		Enters log view and download screens.
			Turns stirrer on and off.

## NOTES

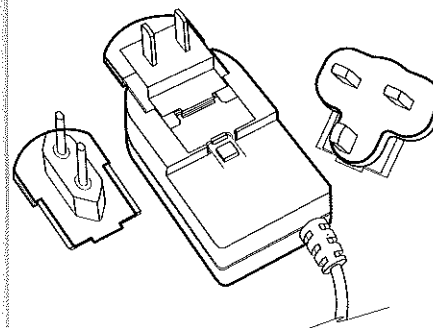
## Chapter IV Preparation

### Installing Power Adapter

The universal power supply that is included with your benchtop meter is the **ONLY** power supply source recommended for use with this unit. Use of any other power supply will void your meter warranty.

The external electrical power supply is rated to be operated at 100-240 VAC, 0.5 A, 50/60 Hz.

Based on your power source, select one of the four plugs provided – 110 V, 220 V, 240V – and slide it into the grooves on the adapter. A click will be heard when the plug is properly in place.




Connect the output plug of the power supply to the power receptacle on the benchtop meter. Refer to the diagram in section **Connecting the Electrodes**.

## Installing Batteries

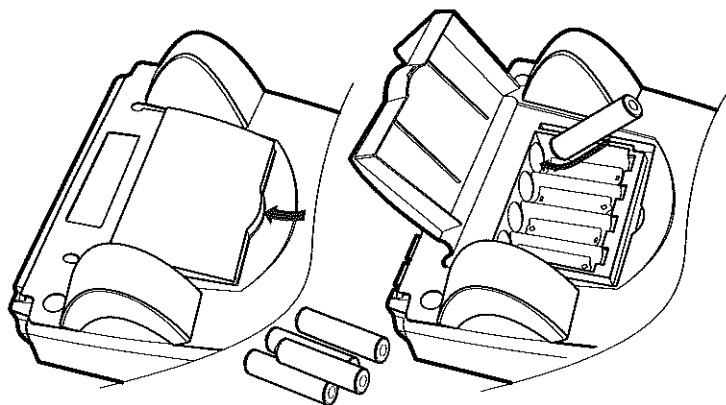The symphony meters use four AA Alkaline batteries. Do not use lithium or rechargeable batteries. Improper installation of non-alkaline batteries could create a hazard.

**Note: For benchtop meters** — installation of batteries is not required if the unit will always be connected to a power source via the universal power supply. ▲


1. Confirm meter is off.
2. Gently place meter upside down on a clean, lint-free cloth to prevent scratching of LCD.
3. Remove the battery case cover.
4. Insert new batteries with the "+" pole oriented as depicted in the battery compartment housing.
5. Replace the cover.



6. Stored data, calibrations and methods will remain in the meter's nonvolatile memory when batteries are being replaced, however, date and time may need to be reset.

**Note: For portable meters** — Batteries are supplied and installed from the factory. ▲

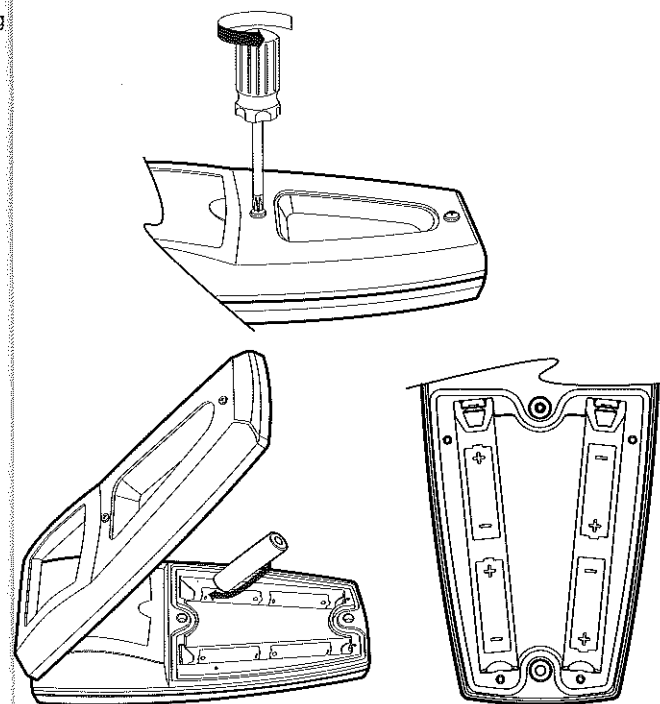To access the battery compartment in portable meters, loosen the two screws in the center back of the meter. Note these are captive screws that cannot be completely removed.







## Connecting the Electrodes

Follow the diagrams below to correctly connect electrodes and probes to the meter. The multi-parameter meter is depicted; dual and single parameter meters will have fewer connections.
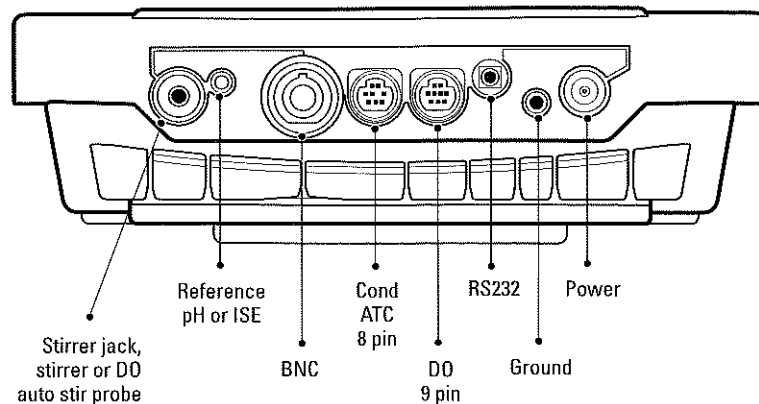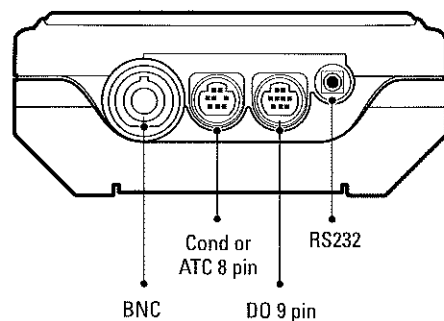
### Benchtop Meter – Electrode Connections







### Portable Meter – Electrode Connections



### Some connections serve multiple uses, for example:

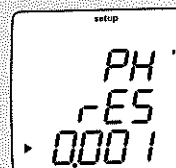
- Use the BNC input to connect pH, ISE and ORP combination and sensing electrodes with a waterproof BNC connector.
- Benchtop meters have a pin tip input for a separate reference electrode. These require an appropriate BNC sensing electrode for measurement.
- The 970899WP dissolved oxygen probe can be used with the BNC connector.
- Use the waterproof 8 pin MiniDIN connector for conductivity probes.
- The waterproof 8 pin MiniDIN connector is also used for automatic temperature compensation (ATC) probes.
- The DO AUTO-STIR™ probe uses the waterproof 9 pin MiniDIN connector and the stirrer jack for the smaller connector.

## Turning On the Instrument



With the batteries installed in the portable meters, or either the batteries installed or line power attached to the benchtop meter, press the  to turn on the instrument. A quick press of the  after the meter has powered up will toggle the backlight on and off. When the benchtop meter is drawing line power, the backlight will stay on until turned off with the . To turn off the meter press and hold the  for 3 seconds.

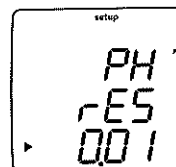
## Chapter V Menu Setup




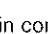


### Navigation Tips



- The  enters the setup mode.

- Main menu items are shown in the top line of the LCD. To scroll through the menu, press  / . Text that will appear in the LCD is represented in the **Display** column of the tables on following pages.



- Use  to select top, middle or bottom line.
- Use  /  in combination with  to edit the values in selected line.
- Use  to accept changes and return the arrow icon to the top line.
- Use  to save changes and return to the measurement mode.

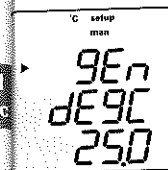
## General Menu Settings

The following table describes general instrument setup.

Description	Display	Range of Values	Default Value	Method Specific
Manual Temperature Compensation Setting	gEn dEgC 25.0	-5.0 - 105	25.0	Yes
Stirrer Speed Setting (Benchtop only)	gEn Stlr 4	OFF, 1, 2, 3, 4, 5, 6, 7	4	Yes
Auto-Shut Off Selection	gEn AUtO On	OFF, On	On	No

- **Manual Temperature** setting controls temperature compensation when no temperature sensors are attached to the instrument.
- **Stirrer Speed Setting** (benchtop meters only) sets the speed from 1 through 7 with speed 1 being the slowest and 7 being the fastest. A speed of 3 or 4 is recommended for most applications.
- **Auto-Shut Off** controls whether the instrument will automatically turn off the instrument after 20 minutes elapse without a key press. To disable this feature, select OFF and the instrument will continue to operate.

To access these settings:

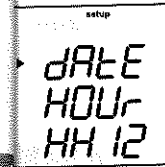














1. In the measurement mode press .
2. Press / to scroll through the setup menu until **gEn** is displayed on the top line.
3. Press to accept selection and move the arrow icon to the middle line.
4. Press / to scroll through:
  - a. **dEgC** = Manual Temperature Setting
  - b. **Stlr** = Enable and Set Stirrer Speed
  - c. **AUto** = Auto-Shut Off enable/disable
5. Press to accept the selection and move the arrow icon to the bottom line.
  - a. Press / to identify and to edit next value.
  - b. Press to accept the selection and move the arrow icon to the top line.
6. Press to return to the measurement mode.

## Time and Date Setup

The following table describes in greater detail the abbreviations shown on the display screen.

Description	Display	Range of Values	Default Value	Method Specific
Time: Hour Setting	dAtE HOUr HH12	00 - 23 Sets the hour of the current time in 24 hour format	12	No
Time: Minute Setting	dAtE tInE mm12	00 - 59 Sets the minutes for the current time	00	No
Date Format	dAtE tYPE dmY	mdY, dmY Month, day, year or day, month, year display format	mdY	No
Date: Year Setting	dAtE YEAr 2004	00 - 99 Sets the year for the current date from 2000 to 2099	04	No
Date: Month Setting	dAtE dAtE mm01	01 - 12 Month setting 01 (January) through 12 (December)	01	No
Date: Day of the Month Setting	dAtE dAY dd01	01 - 31 Day of the month	01	No



1. In the measurement mode press .
2. Press  /  to scroll through the setup menu until **dAtE** is displayed on the top line.
3. Press  to confirm setting and move the arrow icon to the middle line.
4. Press  /  to scroll through:
  - a. **HOUr** = Current Hour
  - b. **tInE** = Current Minute
  - c. **tYPE** = Format the date stamp as either "mdY" (month, day, year) or "dmY" (day, month, year)
  - d. **dAtE** = Current Month (Use numeric equivalents for months of the year)
  - e. **dAY** = Current Day of the Month
  - f. **YEAr** = Current Year
5. Press  to accept selection and move the arrow icon to the bottom line.
6. Press  /  and  to edit the selected value.
7. Press  to accept the selection and move the arrow icon to the top line.
8. Repeat steps 3-7 to edit time and date settings as needed.
9. Press  to return to the measurement mode.

## Continuous, Timed, or AUTO-READ™ Measurement Selection

The following table describes in greater detail the abbreviations shown on the display screen.

Description	Display	Range of Values	Default Value	Method Specific
Continuous, Timed, or AUTO-READ Measurement Selection	rEAd tYPE COnt	<b>COnt, AUtO, tImE</b> Timed, continuous, or AUTO-READ; automatically prints and logs in AUTO-READ and timed modes	<b>AUto</b>	Yes
Timed Reading Setting	rEAd tImE 00:00	<b>00:05 - 99:61</b> Timed readings in minutes and seconds	<b>01:00</b>	Yes

To change settings, simply use and to move between options. Pressing will confirm the change of a value and end the flashing display.

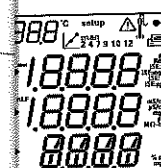
The symPHony meters allow you to select your preferred timing for measurement intervals.



- In continuous mode the instrument is constantly making measurements and updating the screen. To print and log the measurement in this mode press .
- In timed mode the instrument is constantly making measurements, updating the display, and will automatically log and print data at a pre-selected time interval.
  - In the setup menu enter the minutes and seconds between timed readings by pressing / to change the value of the flashing digit. Press to change the digit to be scrolled.
  - The minimum time interval is 5 seconds. Acceptable time ranges are: **00:05 - 99:59**

- In AUTO-READ mode the instrument starts a measurement whenever is pressed. Once the reading becomes stable, the display freezes while data is automatically logged and printed.
- The AUTO-READ mode will also automatically operate the **Stirrer Control**. It will stop stirring when the reading becomes stable. To set the stirrer speed, press . Use / to select:

9Er  
SE Ir  
OFF-7



**AR** AUTO-READ

## Selecting Measurement Mode

In the measurement mode the arrow icon on the left side of the display screen indicates which mode is selected. Use  $\Delta$  /  $\nabla$  to scroll through the various modes associated with a selected line. Use  $\equiv$  to move the arrow icon to the next line, and then toggle  $\Delta$  /  $\nabla$  to scroll through the modes associated with the selected line.

**pH**  
**ISE**  
 mV  
 Rel mV  
 ISE  
 Off

**Cond**  
 $\mu\text{S/cm}$  or  $\text{mS/cm}$  for conductivity  
 mg/L for TDS  
 ppt for Salinity  
 $\text{M}\Omega\text{-cm}$  for resistivity  
 Off

**DO**  
 %Sat for DO percent saturation  
 mg/L for DO concentration  
 Barometric Pressure  
 Off

For setup information related to specific techniques, refer to **Chapter VI** through **Chapter IX** of this manual.

## Setting Calibration

The  $\mathcal{L}$  icon indicates calibration mode or calibration setup.

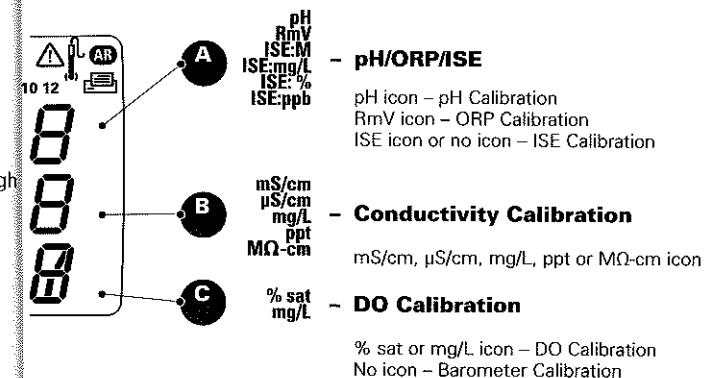
Before beginning a calibration, prepare the electrode or probe according to the instructions received with the electrode or probe. The calibration standards, buffers, and calibration sleeve should also be prepared prior to the calibration. After preparing the electrode or probe per manufacturer's instructions, connect the electrode or probe to the meter. Basic information on proper calibration for each measurement technique is found in the appropriate section of this manual.

### General Navigation for Selecting Calibration

In measurement mode press  $\equiv$  until the arrow icon on the left of the screen is pointed toward the measurement mode to be calibrated:

#### pH, ORP, ISE • Conductivity • DO

Press  $\Delta$  /  $\nabla$  until the appropriate icon is lit for the measurement mode to be calibrated.



Then press  $\mathcal{L}$  to begin the selected calibration.

## NOTES

## Chapter VI pH Technique










### pH Setup

The following table describes pH setup.

Description	Display	Range of Values	Default Value	Method Specific
pH Resolution	PH rES 0.01	<b>0.1, 0.01, 0.001</b>	<b>0.01</b>	Yes
pH Buffer Set	PH bUF USA	<b>USA</b> (USA Buffers) 1.68, 4.01, 7.00, 10.01, 12.46 <b>EUro</b> (Euro Buffers) 1.68, 4.01, 6.86, 9.18	<b>USA</b>	Yes




1. In the measurement mode press
2. Press / to scroll through the setup menu until **PH** is displayed on the top line.
3. Press to confirm setting and move the arrow icon to the middle line.
4. Press / to select **rES** for resolution.
5. Press to select option and move the arrow icon to the bottom line.
  - a. Use / to set the desired resolution.


6. Press  to accept the selection and move the arrow icon to the top line.
7. Press  to move to middle line and press  /  to select **BUF** for Auto-Buffer-Recognition setting and press  again to move to the bottom line.
8. Press  /  to select either **USA** or **EUR-D**; press  to accept the selection.
9. Press  to return to measurement mode.

## Calibration

Prepare the electrode for use according to the instructions received with the electrode.




In setup mode select the buffer set of NIST (**USA**) or DIN (**EUR-D**) being used for the proper buffer recognition to occur, if not set previously.



Press  until the arrow icon is pointing to the pH measurement line.

Press .





Rinse the electrode and ATC probe and place into the buffer.

Wait for the **pH** to stop flashing.

- a. Auto-Buffer-Recognition – When the **pH** stops flashing the meter will display the temperature-corrected pH value for the buffer.
- b. Manual Calibration – When the **pH** stops flashing the meter will display the actual buffer value read by the pH electrode. Use the  and  /  to change the pH value to the temperature-corrected pH value for the buffer.








Once the correct buffer value is shown on the meter display, press  to proceed to the next calibration point and repeat steps 5 through 6 or press  to save the calibration.

The slope will be displayed before the meter returns to the measurement mode. **SLP** is displayed in the lower field and the actual electrode slope, in percent, is displayed in the main field.

- a. On a 1 point calibration use the  and  /  to edit the slope then press  to return to the measurement mode.
- b. On the 2 or more point calibration the meter will automatically go to the measurement mode 2 seconds after display of **SLP**.



## pH Measurement

1. Rinse the electrode in deionized water. Shake off any excess water and blot dry with lint-free tissue.
2. Place the electrode in your sample.
  - a. If you are in continuous measurement mode, the instrument will start measuring immediately. If you are using the benchtop meter and the stirrer control is disabled, then pressing  will start the stirrer. The pH will fluctuate until the reading is stable. Once the reading is stable you can log and print the measurement by pressing . If the stirrer is on, press  to turn off the stirrer before removing the stirrer from the sample.
  - b. If you are in AUTO-READ™ mode then press  to start the measurement. Once the reading is stable the meter will automatically log and print the results and freeze the display. If the stirrer is enabled then the stirrer will turn on when  is pressed and turn off once the reading has stabilized.
  - c. If you are in timed measurement mode then the meter will start making measurements at the frequency you have specified during setup. It will automatically log and print each measurement. If you are using the benchtop meter and the stirrer control is enabled, then pressing  will start the stirrer. Pressing  again will turn off the stirrer.
3. Remove the electrode from the sample and rinse with deionized water, shake off excess water, blot dry, then place it in your next sample and repeat step 2.
4. Once you are finished measuring samples, rinse the electrode with deionized water and blot dry. Then consult your electrode manual for proper electrode storage instructions.

## Chapter VII













# Dissolved Oxygen Technique

## Dissolved Oxygen Setup

The following table describes dissolved oxygen setup.

Description	Display	Range of Values	Default Value	Method Specific
Dissolved Oxygen Saturation Resolution	dO rES 0.1	<b>1, 0.1</b> DO % saturation resolution	<b>0.1</b>	Yes
Dissolved Oxygen Concentration Resolution	dO rES 0.01	<b>0.1, 0.01</b> DO concentration "mg/L" resolution	<b>0.01</b>	Yes
Dissolved Oxygen Barometric Pressure Compensation Selection	dO bAr AUtO	<b>AUto, mAn</b> Selects internal barometer or manual pressure	<b>AUto</b>	Yes
Dissolved Oxygen Manual Barometric Pressure Setting	dO PrES 760.0	<b>450.0 - 850.0</b> Manual pressure compensation value	<b>760.0</b>	Yes
Dissolved Oxygen Salinity Correction Selection	dO SAL AUtO	<b>AUto, mAn</b> Selects method of salinity correction; meters with conductivity function only	<b>AUto</b>	Yes
Dissolved Oxygen Manual Salinity Correction Factor	dO SALF 0	<b>0 - 45</b> Manual salinity correction factor	<b>0</b>	Yes
Dissolved Oxygen Calibration Type Selection	dO CALt Air	<b>Air, H2O, mAn, SEt0</b> Air = Water sat. air H2O = Air sat. water mAn = Manual SEt0 = Zero point	<b>Air</b>	Yes

## DO Setup Steps

1. From the measurement mode, press .
2. Press  /  to scroll through the setup menu until **dO** is displayed on the top line.
3. Press  to confirm the setting and move the arrow icon to the middle line.
4. Press  /  to scroll through:
  - a. **rES + % Sat** = % Saturation Resolution
  - b. **rES + mg/L** = Concentration Resolution
  - c. **bAr** = Barometer Type (Auto/Manual)
  - d. **PrES** = Manual Pressure Compensation
  - e. **SAL** = Automatic/Manual Salinity Compensation
  - f. **SALF** = Manual Salinity Correction
  - g. **CALt** = Calibration Type
5. Press  to select the option and move the arrow icon to the bottom line.
6. Press  /  and  to enter the value. The previous table identifies what you will see on the display and the range of values which may be entered.
7. After entering a value, press  to complete programming of that option and move the arrow icon to the top line. Repeat steps 3 through 6 for all options.
8. Press  to return to measurement mode.

## Calibration


For calibration, the probe must be prepared and polarized.



The DO probe is continuously polarized when connected to the meter. When first connected, or if more than 60 minutes has elapsed with the probe disconnected, re-connect the probe and allow 30 to 60 minutes for polarization. If the probe readings are stable, probe disconnections of less than one hour will require 5 to 25 minutes for polarization.


**Note:** The symphony meters will supply a polarization current to the DO probe even with the power off. To maximize battery life in the meter, unplug the DO probe if it will not be used for an extended period. ▲



Zeroing the probe – a DO probe can generate a 0.02 to 0.05 mg/L positive error in an oxygen free (anoxic) solution. If this error is unacceptable, then zero the probe when using a new sensing membrane, using fresh filling solution, or when measuring dissolved oxygen below 1 mg/L or 10% saturation.


Air calibration should be done prior to the zero calibration.



In the measurement mode press .

Press  /  until **dO** is shown on the top line.



Press  to confirm the setting and move the arrow to the middle line.

Press  /  to scroll to **CALt**.




Press  to select setting and move arrow icon to the bottom line.

Press  /  to select one of the following calibration modes.




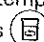
- a. **Air** denotes calibration performed in water saturated air using the air calibration sleeve. This is the simplest and most accurate calibration method.



It is also the meter's default setting. If the calibration selection in  has not been changed, then pressing  will automatically perform an air calibration. Due to the inherent differences between water-saturated air and air-saturated water, upon stability, the air standard is set to 102.3% saturation.

- i. The highest possible accuracy is reached when calibration temperature is the same as the measuring temperature.
  - ii. Moisten the sponge or absorbent cloth in the calibration sleeve with distilled water and insert the probe into the sleeve, but without touching the water saturated material. For BOD measurements, this calibration can be performed in a BOD bottle.
  - iii. For oxygen levels below 1 mg/L, a zero calibration point is often required.
- b. **H<sub>2</sub>O** denotes calibration performed using a water sample that is 100% saturated with air. It is the least commonly used calibration method.
- c. **mAn** denotes manual calibration using a water sample with a known concentration of dissolved oxygen. It can be used to calibrate the sensor to the value achieved by a Winkler titration.
- i. A Winkler calibration involves performing a manual Winkler titration and then using that sample as a standard. The titration oxygen level result is entered in a Winkler calibration as the DO standard value. This correlates the meter input to the Winkler titration. Note, this method is inherently less accurate due to the possibility of titration errors introduced when the calibration is set to the titration test results.
- d. **SELO** denotes a zero calibration, which is used for very low level DO measurements. It is not generally required unless you are making measurements below 5% Saturation or 0.5 mg/L.



7. Press  to accept selection and return the arrow icon to the top line.
8. Press  to return to measurement mode.
9. The probe and calibration standard (water-saturated air or air-saturated water) should be allowed to reach equilibrium before the system is calibrated.
  - a. Press .
  - b. Wait for reading to stabilize.
  - c. Meter will display 102.3% and return to measurement mode.

## DO Measurement

1. Rinse the DO probe, as well as stirrer if being used, in deionized water. Blot dry with a lint-free tissue.
  2. Place the DO probe in your sample.
    - a. If you are in continuous measurement mode, the instrument will start measuring immediately. If you are using the benchtop meter and the stirrer control is enabled, then pressing  will start the stirrer.
- The  will flash until the reading is stable. Once the reading is stable you can log and print the measurement by pressing . If you are recording the data into a notebook, record the DO as well as the temperature at which the DO reading is taken. If the stirrer is in motion, press  to turn off the stirrer before removing the stirrer from the sample.

- b. If you are in AUTO-READ™ mode then press  to start the measurement. Once the reading is stable the meter will automatically log and print the results and freeze the display. If the stirrer is enabled, then the stirrer will turn on when  is pressed and turn off once the reading has stabilized.

When using the AUTO-STIR™ BOD probe, pressing the button on the probe will start the AUTO-READ measurement.





- c. If you are in timed measurement mode then the meter will start making measurements as soon as it goes into measurement mode at the frequency selected in setup. It will automatically log and print each measurement. If you are using the benchtop meter and the stirrer control is enabled then pressing  will start the stirrer. Pressing  again will turn off the stirrer.
3. Remove the DO probe from the sample and rinse with deionized water, then place it in your next sample and repeat step 2.
4. When all samples have been measured, rinse the DO probe with deionized water and blot dry. Consult your DO probe manual for proper electrode storage instructions.



## Chapter VIII

# Conductivity Technique







### Conductivity Setup

Description	Display	Range of Values	Default Value	Method Specific
Conductivity Temperature Compensation Selection	COnD tC Lin	<b>OFF, Lin, nLF</b> Selects temp comp off, linear, or non-linear for natural or ultra-pure water	<b>Lin</b>	Yes
Conductivity Linear Temperature Compensation Coefficient Setting	COnD COEF 2.1	<b>0.0 - 10.0</b> Linear temp compensation coefficient in %/C	<b>2.1</b>	Yes
Conductivity TDS Factor Setting	COnD tdSF 0.49	<b>0.00 - 10.0</b> TDS factor	<b>0.49</b>	Yes
Conductivity Auto-Calibration Default Cell Constant Setting	COnD CELL 0.475	<b>0.001 - 199.0</b> Cell constant used for conductivity AutoCal	<b>0.475</b>	Yes
Conductivity Temperature Reference Selection	COnD trEF 25	<b>15, 20, 25</b> Conductivity reference temperature in °C	<b>25</b>	Yes
Conductivity Cell Type & Manual Ranging Selection	COnD tYPE Std	<b>PLnr, Std, 1, 2, 3, 4, 5, 6, 7</b> Standard, planar cond cell, or sets manual range 1-7	<b>Std</b>	Yes




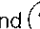

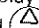
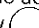

1. In the measurement mode press .
2. Press  /  to scroll through the setup menu until **COnD** is displayed on the top line.
3. Press  to accept the selection and move the arrow icon to the middle line.




4. Press  /  to scroll through the following conductivity settings:

- **TC** = Temperature Compensation Selection. The user may elect to turn temperature compensation **OFF**, select linear **L** temperature compensation, or a non-linear **nLF** temperature compensation for natural and ultra-pure water.
- **COEF** = Temperature coefficient for **L** (linear) temp compensation expressed in % /C.
- **TDSF** = TDS Factor for Total Dissolved Solids measurement
- **CELL** = Conductivity Cell Constant (nominal cell constant used for Auto-Calibration)
- **REF** = Reference temperature for temperature compensation. Options are **15** degrees C, **20** degrees C, or **25** degrees C.
- **TYPE** = Conductivity cell type (Planar or Conventional)


5. Press  to accept the selection and move the arrow icon to the bottom line.
6. Press  /  and  to edit value.
7. Press  to accept the selection and move the arrow to the top line.
8. Repeat steps 3-7 to edit conductivity settings as needed.
9. Press  to return to measurement mode.




## Conductivity Calibration

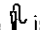



1. Prepare the probe for use according to the probe user's guide or operator's manual.
2. In the measurement mode press  until the arrow icon is pointing to the conductivity measurement line.
3. Press .
4. Rinse the probe and place into the conductivity standard.
5. To perform a manual calibration – The screen will display the cell constant on the bottom line and the conductivity value on the middle line. Use  and  /  to change the cell constant until the displayed conductivity value matches the conductivity of the conductivity standard at the measured temperature.
  - If you do not start changing the cell constant within 5 seconds the meter will automatically proceed to the AutoCal™/DirectCal™ screen.
6. To perform an AutoCal or DirectCal calibration – Wait for the **µS/cm** or **mS/cm** icon to stop flashing.
  - AutoCal – When the **µS/cm** or **mS/cm** icon stops flashing the meter will display the temperature corrected value of the conductivity standard.
  - DirectCal – When the **µS/cm** or **mS/cm** icon stops flashing the meter will display the actual conductivity value read by the conductivity probe. Use  and  /  to change the conductivity value to the actual value of the conductivity standard at the measured temperature.

7. Press  to proceed to the next calibration point and repeat steps 4 through 6 or press  to save and end the calibration.
8. After the value for the last standard has been entered the cell constant will be displayed **CELL** appears in the lower field and the actual **CELL** constant appears in the main field. The meter will automatically advance to the measurement mode. The  is displayed above the main field.

## Conductivity Measurement

1. Rinse conductivity probe and place into sample. Record conductivity directly from the main meter display when the **mS/cm** or **µS/cm** icon stops flashing indicating a stable value. Temperature is displayed in the top left corner of the display.
- 2a. If you are in continuous measurement mode, the instrument will start measuring immediately. If you are using the benchtop meter and the stirrer control is enabled then pressing  will start the stirrer.

The  will flash until the reading is stable. Once the reading is stable you can log and print the measurement by pressing the . If you are recording the data into a notebook, record the conductivity as well as the temperature at which the conductivity reading is taken. If the stirrer is in motion, press  to turn off the stirrer before removing the stirrer from the sample.

- 2b. If you are in AUTO-READ™ mode then press  to start the measurement. Once the reading is stable the meter will automatically log and print the results and freeze the display. If the stirrer is enabled, then the stirrer will turn on when  is pressed and turn off once the reading has stabilized.
- 2c. If you are in timed measurement mode then the meter will start making measurements as soon as it goes into measurement mode at the frequency selected in setup. It will automatically log and print each measurement. If you are using the benchtop meter and the stirrer control is enabled then pressing  will start the stirrer. Pressing  again will turn off the stirrer.





## Conductivity Technique



3. Remove the conductivity probe from the sample and rinse with deionized water, then place it in your next sample and repeat step 2.
4. When all samples have been measured, rinse the conductivity probe with deionized water and blot dry. Consult your conductivity probe manual for proper storage instructions.

## Chapter IX ISE Technique


### ISE Setup




Description	Display	Range of Values	Default Value	Method Specific
ISE Resolution	ISE rES 1	<b>1, 2, 3</b> ISE resolution in significant digits	<b>1</b>	Yes
ISE Units	ISE Unit PPb	<b>M, mG/L, PEr, PPb, nOnE</b> ISE units displayed	<b>PPb</b>	Yes
ISE Calibration Standard Concentration Range	ISE rAng HlgH	<b>LOW, HlgH</b> Stability criteria used during ISE calibration	<b>HlgH</b>	Yes
ISE Auto-Blank Correction	ISE nLin AUtO	<b>AUto, OFF</b> Enables or disables calibration auto-blank	<b>AUto</b>	Yes


1. In the measurement mode press .
2. Press  /  to scroll through the setup menu until **ISE** is displayed on the top line.
3. Press  to accept the selection and move to the middle line.

4. Press  /  to scroll through the following options:


- **rES** = Resolution
- **Un It** = ISE Measurement units
  - **M** = Molar
  - **mG/L** = mg/L
  - **PEr** = Percent
  - **PPb** = Parts per billion
  - **nOnE** = No units
- **rAnG** = ISE calibration range
  - **High** - for most measurements
  - **LOw** - for low level measurements that need more stabilization time
- **nL In** = Non-linear Blank Correction
  - **OFF**
  - **AUtO**

5. Press  to accept the selection and move the arrow icon to the bottom line.

6. Press  /  and  to edit the selected value.

7. Press  to accept the selection and return the arrow icon to the top line.

8. Repeat steps 3 through 7 to edit ISE settings as needed.

9. Press  to return to measurement mode.

## Preparation of Standards

The standards should be prepared in the same ISE units as required by the sample results. It is preferable to use serial dilutions with volumetric glassware to obtain the different concentration levels.








**Note:** Any reagents such as ionic strength adjustors, should be added to samples and standards as specified in the electrode user's guide or instruction manual. ▲

The calibration points should bracket the expected concentration range of the samples and there should be a tenfold change in concentration (i.e. 1 ppm and 10 ppm or 10 ppm and 100 ppm).

Fresh aliquots of standard should be used at each calibration.







## ISE Calibration


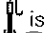

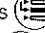

1. Prepare the electrode and standards for use according to the electrode user's guide or instruction manual.
2. Press  until the arrow icon is pointing to the **ISE** measurement line.
3. Press .
4. Rinse the electrode, shake any excess deionized water off and blot dry, and place into the least concentrated standard.
5. Wait for **ISE** to stop flashing. Press  /  and  to change the value of the standard.
6. Press  to proceed to the next lowest calibration standard and repeat steps 4-5, working from lowest concentration to highest, or press  to save the calibration.
7. The slope will be displayed before the meter returns to the measurement mode. **SLP** is displayed in the lower field and the actual electrode slope, in **mV**, is displayed in the main field.








**Note:** For a 2 or more point calibration the meter will automatically proceed to measurement mode after 3 seconds. For a 1 point calibration the meter will allow you to edit the slope then proceed to measurement mode. ▲

**Note:** If editing a negative number for slope: ▲

- a. Press  until no digit is blinking and the arrow icon is blinking.
- b. Press  /  to change the sign of the slope.
- c. Press  to save.

## ISE Measurement

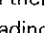After calibrating the electrode the meter is ready to start taking measurements. Make sure that the instrument is in measurement mode ( is lit) and that the **ISE** is lit. If the  is not lit press  to return to measurement mode. If the **ISE** is not lit, press  until the arrow on the left of the display is pointing at the top line. Then press  until the **ISE** lights. You are now ready to start measuring.

1. Rinse the electrode in deionized water. Shake off any excess water and blot dry with lint-free tissue.
2. Place the electrode in your sample.
  - a. If you are in continuous measurement mode, the instrument will start measuring immediately. If you are using the benchtop meter and the stirrer control is disabled, then pressing  will start the stirrer. The **ISE** will flash until the reading is stable. Once the reading is stable you can log and print the measurement by pressing . If the stirrer is on, press  to turn off the stirrer before removing the stirrer from the sample.
  - b. If you are in AUTO-READ™ mode then press  to start the measurement. Once the reading is stable the meter will automatically log and print the results and freeze the display. If the stirrer is enabled then the stirrer will turn on when  is pressed and turn off once the reading has stabilized.
  - c. If you are in timed measurement mode then the meter will start making measurements at the frequency you have specified during setup. It will automatically log and print each measurement. If you are using the benchtop meter and the stirrer control is enabled, then pressing  will start the stirrer. Pressing  again will turn off the stirrer.

## ISE Technique

3. Remove the electrode from the sample and rinse with deionized water, shake or blot dry, then place it in your next sample and repeat step 2.
4. Once you are finished measuring samples, rinse the electrode with deionized water and blot dry. Then consult your electrode manual for proper electrode storage instructions.

## Chapter X Declaration of Conformity

**Address:** 166 Cummings Center  
Beverly, MA 01915  
USA

We declare that the following products described below conform to the Directive and Standard listed below:

**Product(s):** Meters for Measuring pH, Conductivity, Dissolved Oxygen, and/or ISE  
Benchtop models are rated 100-240 VAC, 50/60 Hz, 0.5 A  
Handheld models use 4 non-rechargeable AA Batteries

### **Benchtop**

SB70P benchtop  
SB70D benchtop  
SB70C benchtop  
SB80PI benchtop  
SB80PD benchtop  
SB80PC benchtop  
SB90M5 benchtop

### **Portable**

SP70P portable  
SP70D portable  
SP70C portable  
SP80PI portable  
SP80PD portable  
SP80PC portable  
SP90M5 portable

**Equipment Class:** Measurement, control and laboratory  
Benchtop models are EMC Class A  
Portable models are EMC Class D

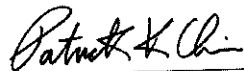
## Declaration of Conformity

### Directive(s) and Standard(s):

- 89/336/EEC – Electromagnetic Compatibility (EMC Directive)
  - EN 61326:1997 + A1:1998 + A2:2001 – Electrical equipment for measurement, control, and laboratory use – EMC requirements
- 73/23/EEC – Low Voltage Directive (LVD)
  - EN 61010-1:2001 – Safety requirements for electrical equipment for measurement, control, and laboratory use – general requirements

Manufacturer's Authorized Representative:

Date:



Patrick Chiu  
Sr. QA Engineering, Regulatory Compliance

June 6, 2007

### WEEE Compliance:



This product is required to comply with the European Union's Waste Electrical & Electronic Equipment (WEEE) Directive 2002/96/EC .  
It is marked with the following symbol:

Thermo Scientific has contracted with one or more recycling/disposal companies in each EU Member State and this product should be disposed of or recycled through them. Further information on compliance with these Directives, the recyclers in your country, and information on Thermo Scientific Orion products which may assist the detection of substances subject to the RoHS Directive are available at [www.thermo.com/WEEERoHS](http://www.thermo.com/WEEERoHS).