


ORNL



NEUTRON SCIENCES

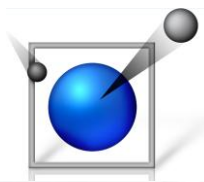

**Fischer Accuspin Micro
17/17R
Microcentrifuge**



Tool for Removing
and Installing
Rotor

14



Ensure the rotor is
balanced with similar
liquid levels across
from each other.

13

Aerosol seals should
be removed if rotor is
used without lid

Power Switch

1



15

Power Cord



Rotation
Display

2

Rotation
Controls

3

RCF/RMP
Selector

4

Time Control

5

Start

6

Stop

7

Open Lid

8

Temperature
Controls

9

Pre-Temp
Mode

10


Rotation
LED

12


Indicator LED
For Temp. Reached

11

SETUP

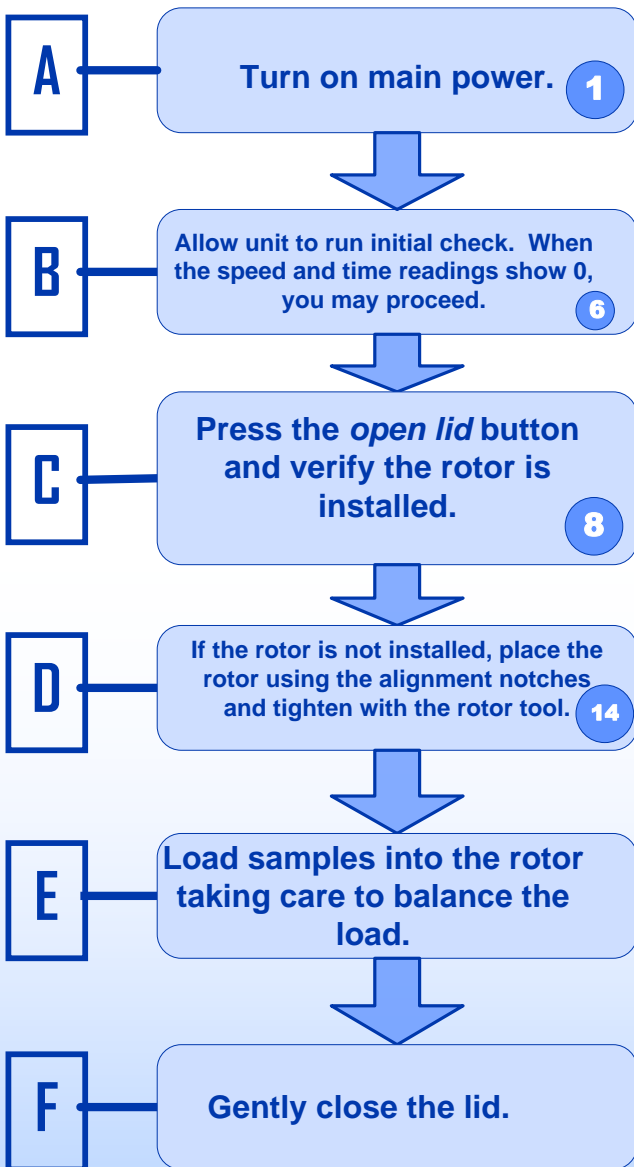

OPERATION

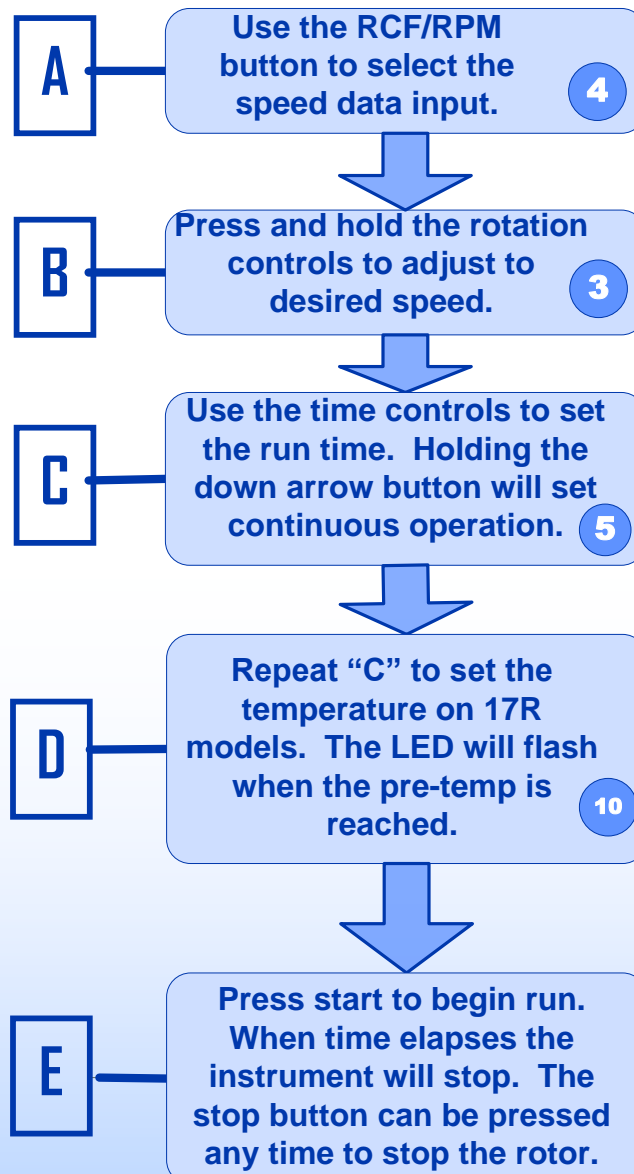

HAZARDS



ELECTRICAL HAZARD

This equipment contains electrical components that may come in contact with liquids or skin.



Slip Hazard

Because liquids are used with this unit, a slip hazard may be created.



PINCH HAZARD

This equipment contains moving parts near the lid/agitator that may cause physical harm.



FIRE HAZARD

Improper equipment use may present a fire hazard.



CHEMICAL HAZARD

Chemical spills or aerosolization may occur.

PREVENTION



EYE PROTECTION

Wear safety glasses while operating.



PROTECTIVE GLOVES

Wear appropriate gloves for heat and chemical usage.



NON-SLIP, CLOSED SHOES

Wear non-slip, closed shoes to avoid spills.



READ MANUAL

Become familiar with this guide before operating.

For additional help contact the Lab
Space Manager