


# INSTRUMENT

BEAM LINE

# 4A


SPALLATION NEUTRON SOURCE

# Fact Sheet

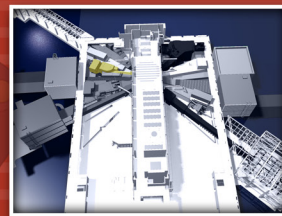

## MAGNETISM REFLECTOMETER

The Magnetism Reflectometer is dedicated to reflectometry studies of magnetic thin films, superlattices, and nanoscale structures. The combination of the high-power SNS and the use of advanced neutron optics allows for off-specular scattering studies of in-plane magnetic and nonmagnetic structures. Wide-angle diffraction geometry is available for experiments on thin films and multilayers. The availability of polarized neutrons and polarization analysis suggests that this instrument can also be used for specific studies of nonmagnetic thin-film samples. Examples of the latter include contrast variation, incoherent background reduction, and phase determination for direct inversion of reflectivity data into real-space scattering-length density profiles.

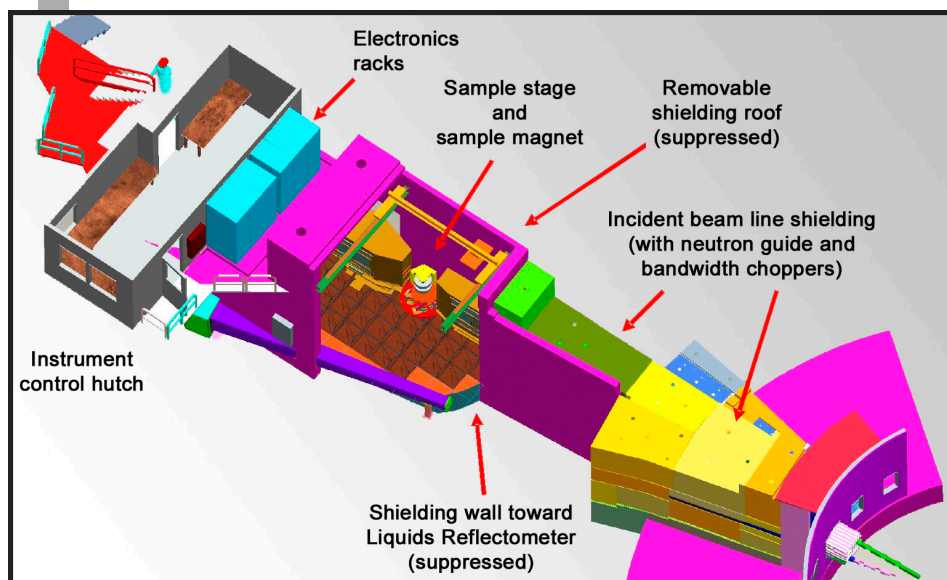

## SPECIFICATIONS

Source-to-sample distance	18.703 m
Sample-to-detector distance	2.540 m
Detector size	21 x 18 cm <sup>2</sup>
Detector resolution	1.5 mm
Moderator	Coupled supercritical hydrogen
Bandwidth	$\Delta\lambda = 2.74 \text{ \AA}$
Wavelength range	$1.8 \text{ \AA} < \lambda < 14.0 \text{ \AA}$
Q range	$0 \text{ \AA}^{-1} < Q < 0.4 \text{ \AA}^{-1}$
Magnetic field max	1.2 T with a gap of 50 mm
Polarization	98.5% and 3 T with a gap of 10 mm
T range	5 - 750 K
Minimum reflectivity	$10^{-8}$

Status: Operational

## APPLICATIONS

The Magnetism Reflectometer is applicable primarily to studies with thin magnetic films, an increasingly important area of solid-state physics. Experiments could also benefit engineering, metallurgy, or biological problems. Instrument capabilities allow, for example, studies of magnetic recording media and magnetic sensors, as well as depth-dependent studies of structural/magnetic nanoparticles or domains. The instrument's unique capabilities provide for multilength-scale experiments, and it has sufficient beam intensity for detailed structural/magnetic phase-diagram determinations.

## FOR MORE INFORMATION, CONTACT

Instrument Scientist: Valeria Lauter, [lauterv@ornl.gov](mailto:lauterv@ornl.gov), 865.387.5389

Scientific Associate: Richard J. Goyette Jr., [goyetterj@ornl.gov](mailto:goyetterj@ornl.gov), 865.241.9991

Scientific Associate: Haile A. Ambaye, [ambayeh@ornl.gov](mailto:ambayeh@ornl.gov), 865.951.8383

<http://neutrons.ornl.gov/instruments/SNS/MR>



May 2011