


INSTRUMENT 4B

BEAM LINE


SPALLATION NEUTRON SOURCE

LIQUIDS REFLECTOMETER

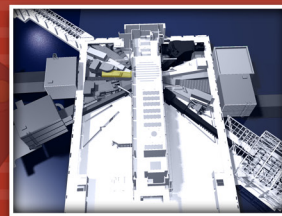

The Liquids Reflectometer features a horizontal sample geometry and thus can accommodate air/liquid surfaces in addition to air/solid and liquid/solid interfaces.

Active vibration isolation minimizes capillary-wave production by the external environment. Data rates and Q range covered at a single scattering angle setting will be sufficiently high to permit “real-time” kinetic studies on many systems. Time-resolved experiments include investigations of chemical kinetics, solid-state reactions, phase transitions, and chemical reactions in general.

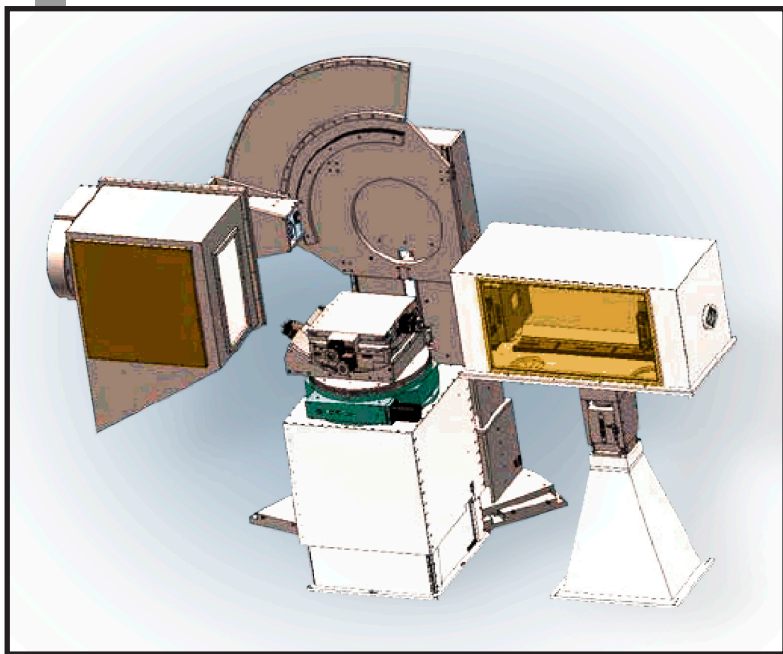

Liquids Reflectometer goniostat.

SPECIFICATIONS

Source-to-sample distance	13.6 m
Sample-to-detector distance	1.5 m
Detector size	20 x 20 cm ²
Detector resolution	1.3 x 1.3 mm ²
Moderator	Coupled supercritical hydrogen
Bandwidth	$\Delta\lambda = 3.5 \text{ \AA}$
Wavelength range	$2.5 \text{ \AA} < \lambda < 17.5 \text{ \AA}$
Q range (air/liquid)	$0 \text{ \AA}^{-1} < Q < 0.4 \text{ \AA}^{-1}$
Q range (air/solid)	$0 \text{ \AA}^{-1} < Q < 0.3 \text{ \AA}^{-1}$
Minimum reflectivity	1×10^{-7}

Status: Operational

APPLICATIONS

The Liquids Reflectometer is useful for a wide range of science. Current areas of interest include biomaterials, polymers, and chemistry involving thin layers of surfactants or other materials on the surfaces of liquids, such as cell-membrane analogs. These systems provide a flexible platform to study structure-property relationships at the boundary between hard and soft matter, with applications in biomimetics, bio-sensing, and bio-compatible films; hydrogen storage and fuel cells; and polymers.

FOR MORE INFORMATION, CONTACT

Instrument Scientist: John Ankner, anknerjf@ornl.gov, 865.377.4523

Instrument Scientist: Jim Browning, browningjf@ornl.gov, 865.209.5273

Scientific Associate: Candice Halbert, halbertce@ornl.gov, 865.643.2334

<http://neutrons.ornl.gov/instruments/SNS/LR>



May 2011